



**SRM VALLIAMMAI ENGINEERING COLLEGE**

(An Autonomous Institution)  
SRM Nagar, Kattankulathur – 603 203



**DEPARTMENT OF CHEMISTRY**

**1901209-Environmental  
Chemistry Lab Manual  
(2022-23)**

**ENVIRONMENTAL CHEMISTRY LABORATORY: (Any five experiments to be conducted)**

**OBJECTIVES**

- To determine the dissolved oxygen and chloride content in water.
- To determine calcium and magnesium present in domestic water.
- To estimate iron, sodium and chlorine using various techniques.
- To determine the chemical oxygen demand in industrial effluent.
- To determine the available chlorine in bleaching powder.

**List of Experiments**

1. Determination of total, temporary & permanent hardness of water by EDTA method
2. Determination of DO content of water sample by Winkler's method
3. Determination of chloride content of water sample by argentometric method
4. Estimation of iron content of the water sample using spectrophotometer
5. Determination of COD value of industrial effluents
6. Estimation of sodium by flame photometry
7. Estimation of available chlorine in bleaching powder

**Demo**

- a. Pollution abatement by adsorption techniques
- b. Scintillation Process

**OUTCOMES**

**The student should be able to:**

- Appreciate the basic requirements for potable water.
- Understand the need of dissolved oxygen in water.
- Explore the quantity of bleaching powder to be added in water.
- Analyze the ill effects caused by the industrial effluents.
- Explore new research areas in the treatment of waste water.

**TEXT BOOK**

G.H. Jeffery, J. Bassett, J. Mendham and R.C. Denney, "Vogel's Textbook of Quantitative Chemical Analysis", John Wiley & Sons Inc, 2014.

## CONTENTS

<b>EX.NO.</b>	<b>NAME OF THE EXPERIMENTS</b>	<b>PAGE NO.</b>
1	Determination of chloride content of water sample by argentometric method	2
2	Determination of total, temporary & permanent hardness of water by EDTA method	6
3	Determination of DO content of water sample by Winkler's method	11
4	Estimation of available chlorine in bleaching powder	15
5	Determination of COD value of industrial effluents	19
<b>EXPERIMENT – BEYOND THE SYLLABUS</b>		
6	Estimation of iron content of the water sample using spectrophotometer	21

**Titration 1 : Standardisation of silver nitrate****Std. NaCl vs AgNO<sub>3</sub>**

Sl. No	Volume of std. NaCl (ml)	Burette reading (ml)		Volume of AgNO <sub>3</sub> (V <sub>1</sub> ml)	Concordant value (ml)	Indicator
		Initial	Final			
1.	20	0				K <sub>2</sub> CrO <sub>4</sub>
2.	20	0				

**Calculation of the Strength of silver nitrate**

Volume of std. NaCl (V<sub>2</sub>) = 20 ml

Strength of NaCl (N<sub>2</sub>) = 0.01 N

Volume of AgNO<sub>3</sub> (V<sub>1</sub>) = \_\_\_\_\_ ml

Strength of AgNO<sub>3</sub> (N<sub>1</sub>) = ?

According to volumetric formula

$$V_1 N_1 = V_2 N_2 \text{ i.e. } N_1 = V_2 \times \frac{N_2}{V_1}$$

$$N_1 = 20 \text{ ml} \times 0.01 \text{ N} / V_1$$

Strength of AgNO<sub>3</sub> (N<sub>1</sub>) = \_\_\_\_\_ N

## 1. ESTIMATION OF CHLORIDE CONTENT IN WATER BY ARGENTOMETRIC METHOD (MOHR'S METHOD)

Expt. No.	
-----------	--

Date			
------	--	--	--

### AIM

To estimate the amount of chloride ion present in the water sample by Argentometric method (Mohr's method). You are provided with standard NaCl solution of strength 0.01 N and a link solution of AgNO<sub>3</sub> (approximately 0.01N).


### CHEMICALS REQUIRED

Standard NaCl solution, AgNO<sub>3</sub> solution, Potassium chromate indicator


### PRINCIPLE

Natural water contains chloride ions in the form of NaCl, KCl, CaCl<sub>2</sub>, MgCl<sub>2</sub>. The concentration of chloride ions more than 250 ppm is not desirable for drinking purposes.

This determination is based on precipitation titration. When AgNO<sub>3</sub> solution is added to the water sample, in presence of K<sub>2</sub>CrO<sub>4</sub>, the chlorides present in it are precipitated first as AgCl.

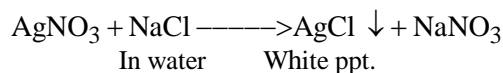

When all the Cl<sup>-</sup> ions is removed, AgNO<sub>3</sub> added from the burette will react with K<sub>2</sub>CrO<sub>4</sub> to give a reddish brown colour due to silver chromate (Ag<sub>2</sub>CrO<sub>4</sub>). This is the end point.

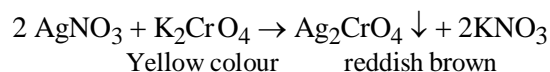

### PROCEDURE

#### Titration 1: Standardisation of AgNO<sub>3</sub>

The burette is washed well with distilled water and rinsed with small amount of AgNO<sub>3</sub> solution. The pipette is washed with distilled water and rinsed with small amount of standard NaCl solution. 20ml of this solution is pipetted out into a clean conical flask. 1ml of 2% K<sub>2</sub>CrO<sub>4</sub> indicator solution is added and titrated against AgNO<sub>3</sub> solution taken in the burette. The end point is the change of colour from yellow to reddish brown. The titration is repeated for concordant values.

**Titration 2 : Estimation of Chloride****Water sample vs Standard AgNO<sub>3</sub>**

Sl. No	Volume of water sample (ml)	Burette reading (ml)		Volume of AgNO <sub>3</sub> (ml)	Concordant value (ml)	Indicator
		Initial	Final			
1.	20	0				K <sub>2</sub> CrO <sub>4</sub>
2.	20	0				

**Calculation of the normality of water sample (Chloride ion)**

Volume of water sample (V<sub>1</sub>) = 20 ml

Strength of water sample (N<sub>1</sub>) = ?

Volume of AgNO<sub>3</sub> (V<sub>2</sub>) = \_\_\_\_ ml

Strength of AgNO<sub>3</sub> (N<sub>2</sub>) = \_\_\_\_ N

According to volumetric formula

$$V_1 N_1 = V_2 N_2$$

$$N_1 = V_2 \times N_2 / 20$$

Strength of water sample (N<sub>1</sub>) = \_\_\_\_ N

**Calculation of amount of chloride**

Amount of chloride present in 1 litre of the given water sample = Eq.wt of chloride ion × Normality of chloride ion

Amount of chloride ion present in 100 ml of the given water sample =  $35.46 \times \text{Normality of chloride ion} \times 100/1000\text{gm}$

**Titration 2: Estimation of Chloride ion**

The given water sample is made up into 100 ml in a standard flask using distilled water. 20ml of this solution is pipetted out into a clean conical flask and 1ml of 2%  $K_2CrO_4$  indicator solution is added. It is then titrated against standard  $AgNO_3$  solution taken in the burette. The end point is the change of colour from yellow to reddish brown. The titration is repeated for concordant values.

**RESULT**

The amount of chloride ion present in 100 ml of the given water sample = \_\_\_\_\_gms.

**Titration: 1 Determination of Total Hardness****Sample hard water vs standard EDTA**

S. No	Volume of sample hard water (ml)	Burette reading (ml)		Volume of EDTA (ml)	Concordant value (ml)	Indicator
		Initial	Final			
1.	20	0				Eriochrome black -T
2.	20	0				

Volume of hard water sample ( $V_1$ ) = 20ml

Strength of hard water sample ( $N_1$ ) = N

Volume of EDTA ( $V_2$ ) = — ml

Strength of EDTA ( $N_2$ ) = 0.01N

According to volumetric formula

$$V_1 N_1 = V_2 N_2$$

$$N_1 = \frac{V_2 \times 0.01}{20}$$

Strength of hard water sample ( $N_1$ ) = — N

Total hardness of the hard water sample =  $N_1 \times 50 \times 1000$

= ppm

VEC - Chemistry



## 2. ESTIMATION OF HARDNESS OF WATER BY EDTA METHOD

Expt. No.	
-----------	--

Date			
------	--	--	--

### AIM

To estimate the amount of total hardness, permanent hardness and temporary hardness of a given hard water sample by EDTA method. You are provided with a standard solution of calcium of strength 0.01 M.

### CHEMICALS REQUIRED


EDTA, Eriochrome black-T, water sample, ammoniacal buffer solution.

### PRINCIPLE

The amount of calcium and magnesium ions present in a given sample of water is estimated by titration with standard EDTA. The hardness is expressed in terms of calcium carbonate equivalence.

### Hardness

It is one of the important characteristics of water. Hardness of water is due to the presence of calcium and magnesium salts. Lather is not produced until the cations,  $\text{Ca}^{2+}$  and  $\text{Mg}^{2+}$  are completely precipitated out in the form of insoluble salts. For example,

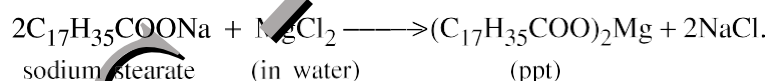

### Temporary hardness

It is due to bicarbonates of calcium and magnesium. It can be removed by boiling. The salts are precipitated as carbonates and hydroxides respectively, which can be removed by filtration leaving only permanent hardness producing salts.

### Permanent hardness

It is due to the sulphates and chlorides of calcium and magnesium.

### Hardness is expressed in terms of $\text{CaCO}_3$ equivalence

The reason for choosing  $\text{CaCO}_3$  as the standard for calculating hardness of water is due to its molecular weight being 100, which makes the calculations easier.

- It is the most insoluble salt, thus can be easily precipitated in water treatment processes.

**Titration: 2 Determination of Permanent Hardness****Boiled water sample vs Standardised EDTA**

S. No	Volume of boiled hard water Sample (ml)	Burette reading (ml)		Volume of EDTA (ml)	Concordant value (ml)	Indicator
		Initial	Final			
1.	20	0				Eriochrome black -T
2.	20	0				

Volume of boiled hard water sample ( $V_1$ ) = 20ml

Strength of boiled hard water sample ( $N_1$ ) = N

Volume of EDTA ( $V_2$ ) = --- ml

Strength of EDTA ( $N_2$ ) = 0.01N

According to volumetric formula

$$V_1 N_1 = V_2 N_2$$

$$N_1 = \frac{V_2 \times 0.01}{20}$$

$$N_1 = \frac{\quad}{20}$$

Strength of boiled hard water sample ( $N_1$ ) = --- N

Permanent hardness of the hard water sample =  $N_1 \times 50 \times 1000$

= ppm


Temporary hardness = Total hardness - Permanent hardness

= .....ppm


### Determination of hardness

The hardness of water can be determined by complexometric titration. EDTA is used as a complexing agent. The  $\text{Ca}^{2+}$  and  $\text{Mg}^{2+}$  present in water are titrated with EDTA using Eriochrome Black T as indicator.

The Chemical structure of EDTA (Ethylene diamine tetra acetic acid) is represented as below:

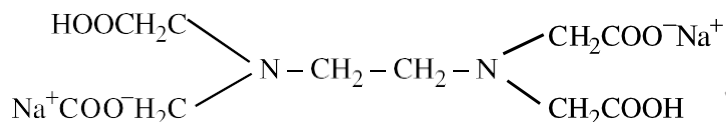


### STRUCTURE OF EBT

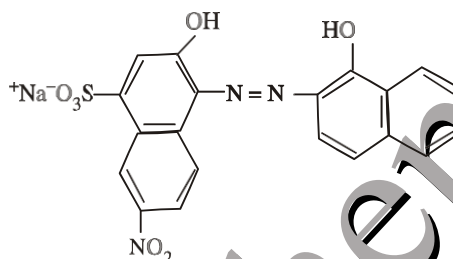

The indicator Eriochrome Black-T (EBT), forms a wine red colored unstable complex with calcium and magnesium ions in hard water at a pH of 8-10.



As this solution is titrated against EDTA, the  $\text{Ca}^{2+}$  and  $\text{Mg}^{2+}$  ions from the indicator complex forms stable metal ion- EDTA complex leaving the indicator free. The free form of indicator is steel blue in colour. The metal-EDTA complex is stable at pH 8-10. This pH range can be maintained by using ammoniacal buffer ( $\text{NH}_4\text{Cl} + \text{NH}_4\text{OH}$ ).

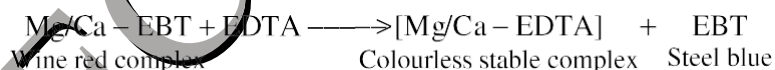

Thus the amount of EDTA used corresponds to the hardness of water.

**PROCEDURE :****Titration : 1 Estimation of total hardness of sample water**

20 ml of the given hard water sample is pipetted out into a clean conical flask. 5ml of ammonical buffer solution and 2 drops of EBT indicator are added. The solution turns wine red in colour. The solution is then titrated against EDTA solution taken in the burette. The change of colour from wine red to steel blue is the end point. The titration is repeated for concordant values. Let the titre value be  $V$  ml. This gives the total hardness of the sample water.

**Titration : 2 Estimation of permanent Hardness**

100 ml of the given sample water is taken in a clean 250 ml beaker and boiled for 10-15 minutes until it is reduced to  $1/5^{\text{th}}$  of the original volume. It is then cooled and filtered. The filtrate is collected in a 100 ml standard flask and made upto the mark using distilled water. 20 ml of this made up solution is pipetted out into a clean conical flask and the titration is repeated in the same way as in titration 1. Let the titre value be  $V$  ml. This gives the permanent hardness of the sample water.

**Estimation of temporary hardness**

Temporary hardness of the water sample is calculated by subtracting permanent hardness from total hardness.

$$\text{Temporary hardness} = \text{Total Hardness} - \text{Permanent hardness.}$$

**RESULT**

1. The total hardness of the given water sample = \_\_\_\_\_ ppm
2. The permanent hardness of the given water sample = \_\_\_\_\_ ppm
3. The temporary hardness of the given water sample = \_\_\_\_\_ ppm.

**Step I : Standardization of sodium thiosulphate****Titration I (Standard dichromate vs Sodium thiosulphate)**

Sl.No	Volume of Potassium dichromate (ml)	Burette reading (ml)		Volume of Sodium thiosulphate (ml)	Indicator
		Initial	Final		
1	20	0			Starch
2	20	0			
3	20	0			

Volume of Potassium dichromate ( $V_1$ ) = 20 ml

Strength of Potassium dichromate ( $N_1$ ) = 0.0125 N

Volume of Sodium thiosulphate ( $V_2$ ) = --- ml

Strength of Sodium thiosulphate ( $N_2$ ) = --- N

According to volumetric formula

$$V_1 N_1 = V_2 N_2$$

$$N_2 = \frac{V_1 N_1}{V_2}$$

Strength of Sodium thiosulphate ( $N_2$ ) = --- N

### 3. ESTIMATION OF DISSOLVED OXYGEN OF BOILER FEED WATER

<b>Expt. No.</b>	
------------------	--

<b>Date</b>			
-------------	--	--	--

#### AIM


To determine the dissolved oxygen in the boiler feed water.

#### CHEMICALS REQUIRED




Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub> (N/40), MnSO<sub>4</sub> solution, KI, starch, conc. H<sub>2</sub>SO<sub>4</sub>

#### PRINCIPLE

Oxygen dissolves in water to the extent of 7 - 9 mgs/lit at a temperature range of 25 - 35°C. The estimation of dissolved oxygen in water is useful in studying water pollution. Water sample is collected carefully avoiding aeration/deaeration in ground stoppered flask. Initially manganous sulphate and alkali-iodide reagents are added and the reaction occur as follows

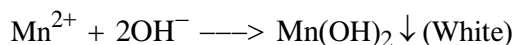

Potassium iodide and the precipitate react with concentrated sulphuric acid liberating iodine and the liberated iodine is titrated against Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub>

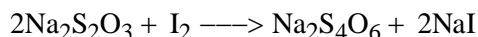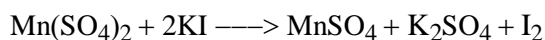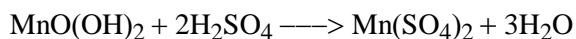

#### PROCEDURE

##### Titration I

##### Standardisation of sodium thiosulphate

The burette is washed and rinsed with sodium thiosulphate solution. Then the burette is filled with the given sodium thiosulphate solution. 20 ml of 0.0125N potassium dichromate solution is pipetted out into a clean conical flask. To this, 5 ml of sulphuric acid and 50 ml of 5% potassium iodide are added. This is titrated against sodium thiosulphate solution. When the solution becomes straw yellow colour, starch indicator is added and then titration is continued.

The end point is disappearance of blue colour and appearance of light green colour. The titration is repeated to get concordant values.

**Step II : Estimation of dissolved oxygen****Titration II (Water sample vs Sodium thiosulphate)**

Sl.No	Volume of Water sample (ml)	Burette reading (ml)		Volume of Sodium thiosulphate (ml)	Indicator
		Initial	Final		
1	100	0			Starch
2	100	0			
3	100	0			

Volume of Sodium thiosulphate ( $V_1$ ) =

Strength of Sodium thiosulphate ( $N_1$ ) =

Volume of water sample ( $V_2$ ) = 100 ml

Strength of water sample ( $N_2$ ) = --- N

According to volumetric formula

$$V_1 N_1 = V_2 N_2$$

$$N_2 = \frac{V_1 N_1}{V_2}$$

Amount of dissolved oxygen in one litre of tap water =

Normality  $\times$  Eq.wt.of  $O_2 \times 1000$  mg.

$$= \text{---} N \times 8 \times 1000$$

$$= \text{---} \text{ mg/lit}$$

## **Titration II**

### **Estimation of dissolved oxygen**

100-150ml of the water sample is taken in the iodine flask, 2ml of manganese sulphate and 2 ml of alkali-iodide are added. The stopper is replaced and the flask is inverted and shaken several times for thorough mixing of the reagents. The flask is left aside for sometime. When half of the precipitate settles down, the stopper is removed and 2 ml of concentrated sulphuric acid is added.

The stopper is replaced and the flask is inverted several times for complete dissolution of the precipitate. 100 ml of the brown coloured solution is pipetted out and titrated against standardized sodium thiosulphate solution. Starch indicator is added when the solution becomes light yellow. The titration is continued until the blue colour disappears. From the titre value, the strength of dissolved oxygen and hence the amount of dissolved oxygen in the water sample is calculated.

### **RESULT**

Amount of dissolved oxygen in water sample = ----- mg/lit



**Titration I: Standardization of sodium thiosulphate****Standard dichromate vs Sodium thiosulphate**

Sl.No.	Volume of Potassium dichromate (ml)	Burette readings (ml)		Volume of sodium thiosulphate (ml)	Concordant Value (ml)
		Initial	Final		
1	20	0			
2	20	0			

Volume of Potassium dichromate ( $V_1$ ) = 20 ml

Strength of Potassium dichromate ( $N_1$ ) = 0.01 N

Volume of Sodium thiosulphate ( $V_2$ ) = --- ml

Strength of Sodium thiosulphate ( $N_2$ ) = --- N

According to volumetric formula,  $V_1N_1 = V_2N_2$

$$N_2 = \frac{V_1N_1}{V_2}$$

Strength of Sodium thiosulphate ( $N_2$ ) = --- N

#### 4. ESTIMATION OF AVAILABLE CHLORINE IN BLEACHING POWDER

<b>Expt. No</b>	
-----------------	--

<b>Date</b>			
-------------	--	--	--

#### AIM:




To determine the Percentage of Available chlorine in the given sample of bleaching powder solution Iodometrically.

#### CHEMICALS REQUIRED:


N/10 Hypo ( $\text{Na}_2\text{S}_2\text{O}_3 \cdot 6\text{H}_2\text{O}$ ), bleaching powder sample solution, Solid KI, dilute  $\text{H}_2\text{SO}_4$ .

#### PRINCIPLE:

The amount of chlorine liberated by the action of dilute acids on bleaching powder is termed as available chlorine and expressed as percentage weight of bleaching powder. The bleaching powder is a mixture of  $(\text{CaOCl}_2 \cdot 4\text{H}_2\text{O})$  ( $\text{CaCl}_2 \cdot \text{Ca}(\text{OH})_2 \cdot \text{H}_2\text{O}$ ) and some free  $\text{Ca}(\text{OH})_2$ . When dilute  $\text{H}_2\text{SO}_4$  reacts with bleaching powder then free chlorine is liberated. The liberated chlorine reacts with potassium iodide to give free iodine. This liberated free iodine is then titrated against N/10 Hypo solution, using freshly prepared starch solution.

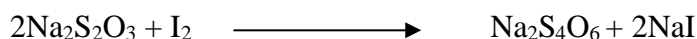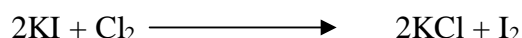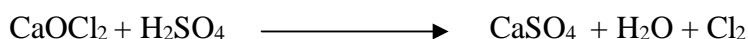

Sodium Tetra Thionate

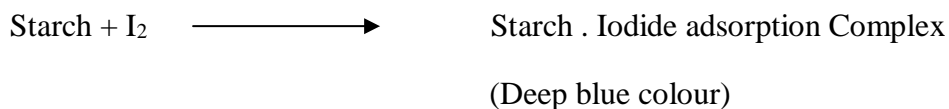

#### PROCEDURE:

Titration 1:

#### Standardization of sodium thiosulphate

The burette is washed and rinsed with sodium thiosulphate solution. Then the burette is filled with the given sodium thiosulphate solution. 20 ml of 0.0125N potassium dichromate solution is pipetted out into a clean conical flask. To this, 5 ml of sulphuric acid and 50 ml of 5% potassium

**Titration 2: Estimation of available chlorine****Bleaching powder vs Sodium thiosulphate**

Sl.No.	Volume of Bleaching Powder sample taken (ml)	Burette readings (ml)		Volume of sodium thiosulphate (ml)	Concordant Value (ml)
		Initial	Final		
1	20	0			
2	20	0			

Volume of Sodium thiosulphate ( $V_1$ ) = ml

Strength of Sodium thiosulphate ( $N_1$ ) = N

Volume of Bleaching Powder solution ( $V_2$ ) = 20 ml

Strength of Bleaching Powder solution ( $N_2$ ) = — N

According to volumetric formula,  $V_1N_1 = V_2N_2$

$$N_2 = \frac{V_1N_1}{V_2}$$

Strength of Bleaching Powder solution ( $N_2$ ) = — N

Amount of Chlorine per Litre of solution = Normality x Equivalent weight

$$= N_2 \times 35.5 \text{ g/lit}$$

Percentage of Available Chlorine =  $N_2 \times 35.5 \times \frac{500}{1000} \times \frac{100}{W}$

= .....%

iodide are added. This is titrated against sodium thiosulphate solution. When the solution becomes straw yellow colour, starch indicator is added and then titration is continued.

The end point is disappearance of blue colour and appearance of light green colour. The titration is repeated to get concordant values.

### **Titration 2:**

#### **Estimation of available chlorine**

Fill the burette with Hypo solution. Pipette out 20 ml of bleaching power solution into 250 ml conical flask. Add 2 ml of KI and 10 ml of dilute  $H_2SO_4$ . Cover the mouth of conical flask with lid and keep it for 2 minutes. Solution becomes brown. The liberated Iodine against Hypo solution till a light yellow colour persists. Add 1 to 2 ml of starch solution. Solutions turns to blue colour. Continue adding Hypo solution till blue colour disappears. This the end point. Repeat the experiment until the concordant readings are obtained.

#### **RESULT:**

The percentage of available chlorine in the given sample of bleaching powder is .....

**Tabulation :**

Sl.No.	Volume of water sample (ml)	Burette readings (ml)		Volume of Ferrous Ammonium Sulphate (ml)	Concordant Value (ml)
		Initial	Final		

**Calculation:**

Quantity of  $\text{Fe}(\text{NH}_4)_2(\text{SO}_4)_2$  added for blank (A) =

Quantity of  $\text{Fe}(\text{NH}_4)_2(\text{SO}_4)_2$  added for sample (B) =

$$\begin{aligned} \text{Dilution Factor} &= \frac{\text{Volume of dilute sample}}{\text{Volume of waste water sample added}} \\ &= \frac{40}{10} \\ &= 4 \end{aligned}$$

$$\text{Chemical Oxygen Demand (COD)} = \frac{(A-B) \times N \times 8 \times 1000 \times \text{Dilution Factor}}{\text{Volume of Sample}}$$

$$\text{COD} = \dots\dots\dots$$

## 5. DETERMINATION OF COD VALUE OF INDUSTRIAL EFFLUENTS

<b>Expt. No</b>	
-----------------	--

<b>Date</b>			
-------------	--	--	--

**AIM:**

To determine the chemical oxygen demand (COD) exerted by the given waste water sample.

**APPARATUS REQUIRED:**

Reflux Apparatus, Burette, Hot plate/ Heating mantel.

**PRINCIPLE:**

The organic matter present in sample gets oxidized completely by  $K_2Cr_2O_7$  in the presence of  $H_2SO_4$  to produce  $CO_2$  and  $H_2O$ . The excess of  $K_2Cr_2O_7$  remaining after the reaction is titrated with  $Fe(NH_4)_2(SO_4)_2$ . The dichromate consumed gives the oxygen required for the oxidation of organic matter.


**REAGENTS:**

1. Standard Potassium Dichromate 0.2N
2. Sulphuric Acid with reagent (Conc.  $Ag_2SO_4 + H_2SO_4$ )
3. Std. Ferrous Ammonium Sulphate (0.1N)
4. Ferroin Indicator
5. Mercuric Sulphate
6. Silver Sulphate


**PROCEDURE:**

1. Take 10 ml of samle.
2. 5 ml of more conc. Dichromate solution are placed in a flask together with glass beads.
3. Add slowly 15 ml of  $H_2SO_4$  containing  $Ag_2SO_4$  and mix thoroughly.
4. Add pinch of mercurous sulphate ( $HgSO_4$ ) and silver sulphate ( $Ag_2SO_4$ )
5. Connect the flask to condenser. Mix the contents thoroughly before heating.  
Improper mixing may result in bumping and the sample may be blown out.
6. Reflux for a minimum period of 2 hours. Cool and wash down the condenser with distilled water.
7. Dilute the sample to make up to 40 ml with distilled and cool.
8. Add 2-3 drops of Ferroin indicator. Mix thoroughly and titrate it against 0.1 N  $Fe(NH_4)_2(SO_4)_2$ . Sharp colour change from blue- green to wine red indicate the end point.
9. Reflux the blank solution in the same manner using distilled water instead of sample.

**RESULT:** The COD for the given water sample is.....

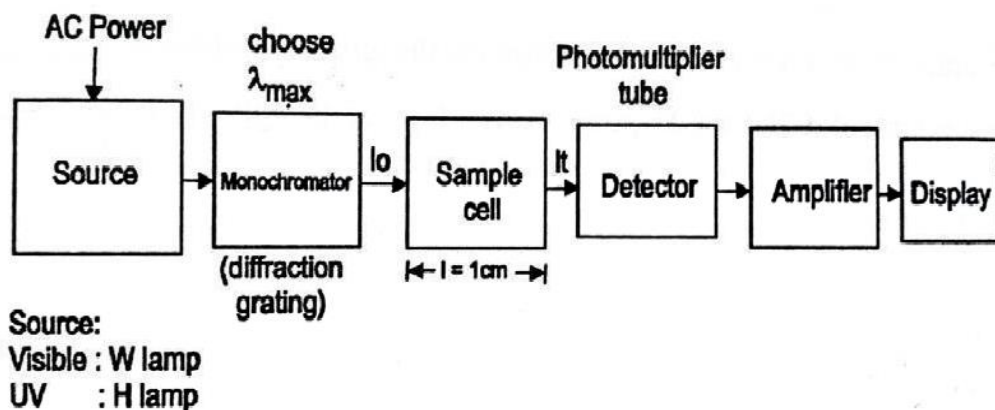

**UV - Vis Spectrophotometer - Block diagram**

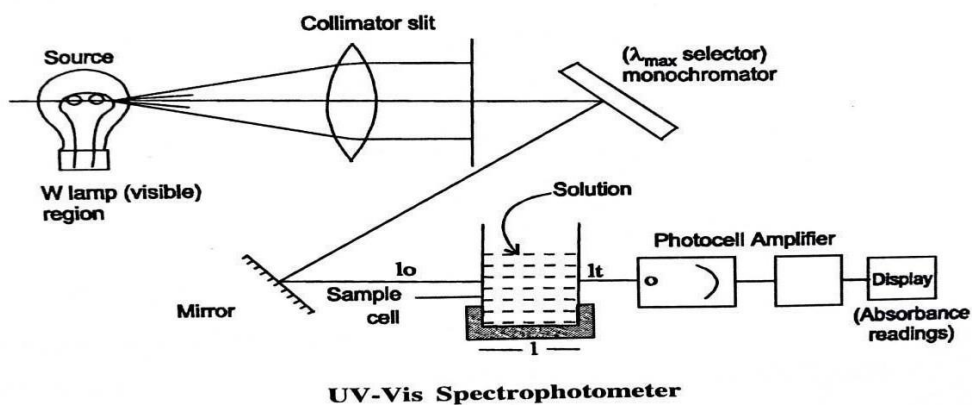

**UV-Vis Spectrophotometer**

## V. DETERMINATION OF IRON IN WATER SAMPLE BY SPECTROPHOTOMETRY (USING THIO CYANATE)

<b>Expt. No</b>	
-----------------	--

<b>Date</b>			
-------------	--	--	--

**AIM:**

To estimate the amount of  $\text{Fe}^{3+}$  ion or  $\text{Fe}^{2+}$  ion present in the given aqueous solution by spectrophotometry. You are provided with a standard solution of ferrous ammonium sulphate containing 0.0838g in 1000ml.

**MATERIAL REQUIRED:**

Spectrophotometer, ferrous ammonium sulphate, nitric acid, potassium or ammonium thiocyanate, glass wares.

**PRINCIPLE:**

When a light passed through a solution, a portion of the light is absorbed. The portion of absorbed light depends on the concentration of the solution,


According to Lamberts Beer Law, if  $I_0$  and  $I_t$  are the intensities of the incident and transmitted light.

$$\text{Log } I_0/I_t = \epsilon cl = A$$

$\epsilon, l$  are constants.  $A$  = absorbance,  $c$  = concentration

Thus by finding out the change in intensity of colour of the light absorbed by a solution, its concentration can be calculated.

In order to get a visible colour,  $\text{Fe}^{2+}$  ion is oxidized to  $\text{Fe}^{3+}$  by nitric acid and treated with potassium thiocyanate to get deep red colour.

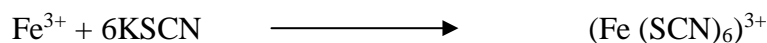

(red colouredthiocyanato complex of iron  
which absorbs in the visible region of 480 – 500 nm)


**PROCEDURE:**

A standard solution may be prepared by dissolving 0.0838 g of pure ferrous ammonium sulphate in distilled water with 1 ml of con.  $\text{HNO}_3$  and making upto 1000 ml (containing 1 ppm of iron / ml). 2,4,6 and 8 ppm iron solutions are prepared by pipetting out 2,4,6 and 8 ml solution in to a measuring cylinder,



**Absorbance of Solution at 480 – 500 nm**

S.No.	Concentration of Solution (mg/l)	Absorbance(nm)
1		
2		
3		
4		
5		

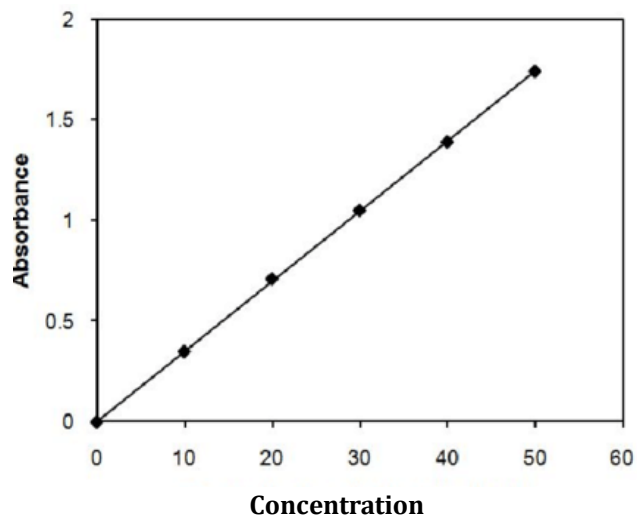**Graph**

1ml of 1:1 HNO<sub>3</sub> and 1 ml of 10% thiocyanate and remaining distilled water are added to make the total volume of 10 ml.

The spectrophotometer is adjusted for  $\lambda = 480\text{nm}$ . The transmittance ( $I/I_0$ ) of the blank solution taken in the cell is adjusted for 100% absorbance. Solutions of various known concentrations are taken one after the other and the respective absorbances are measured. Lastly, the unknown solution is taken in the cell after carrying out the same treatment with nitric acid and thiocyanate.

Absorbance is measured. A graph is drawn connecting absorbance with concentration. A straight line is expected. From the graph and the absorbance of the unknown solution its concentration can be found out.

**RESULT:**

Amount of iron present in the given sample of water is = ----- ppm.